

1. Especifications

Model:	3000835
RPM	60...200
Time:	From 1 minute to 24h or continuous operation (TIME=0)
Power:	110-230VAC 50/60Hz 20W
	12VAC 20W (car cigarette lighter outlet)


2. Instalation

Place on a flat, level surface.

Leave about 30cm free around the equipment to allow manipulation of the vessels.

Use preferably 500ml beakers, low form, Ø110 x 145 mm.

3. Parts

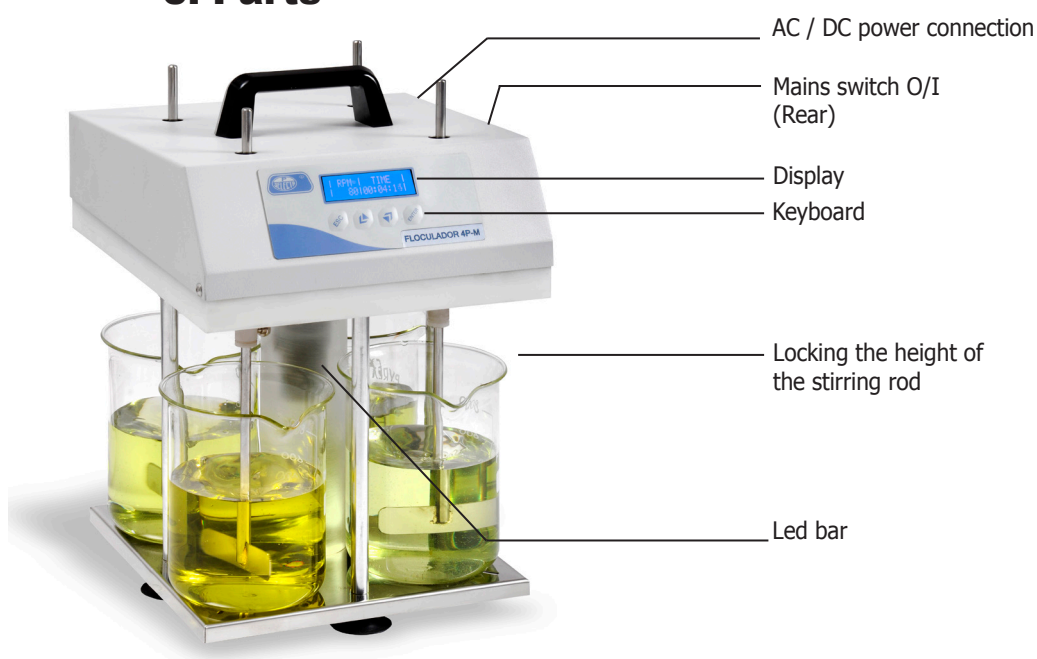

4. Operation

Initial screen:

ESC

Set speed (rpm):



Entry value

ENTER

Set time (hh:mm:ss):
00-00 for continuous operation



Entry value

ENTER

Set language (ESP, ENG, FRA)



Entry value

ENTER

Start operation

ENTER

Stop operation

ESC

(Returns to initial screen)

Alarm sounds if viscosity is too high.