

APPARATUS FOR WATER ANALYSIS



Thermo reactor for determining COD Chemical Oxygen Demand of waste water plants

STANDARD TO C.E.E. AND UNE 77-004.
REFLUX METHOD.

FEATURES

Composed of a Hot block for 6,12 or 20 samples, process programmer and lifting racks. Uniform heating throughout the block to all samples. Automatic temperature and time control during digestion.

C.O.D. unit comprises of:

- 1 metallic hot block.
- 1 Processor for temperature and time.
- 1 Support rack for tubes.
- 1 Support rack for condensers.
- C.O.D. tubes with an aperture of 29/32.
- Condensers for C.O.D.

MODELS - EQUIPMENT DESCRIPTION

	Part No.	No. of places
DQO. - 6	4000638	6
DQO. - 12	4000639	12
DQO. - 20	4000640	20



INDIVIDUAL PARTS LIST AND SPARES

Hot block. This block cannot be used without the RAT controller of temperature and time.




Part No.	No. of places	Ø tube mm	Height / Width / Depth (exterior) cm	Temperature °C	Power W	Weight Kg
4000507	6	42	18 33 28	45 - 450	1500	18
4000508	12	42	18 39 33	45 - 450	2100	25
4000509	20	42	18 44 39	45 - 450	2500	31



System programmer for temperature/time RAT-2. Part No. 4001538

Condenser tube support.

- Part No.
- 4000643 6 places.
 - 4000644 12 places.
 - 4000645 20 places.

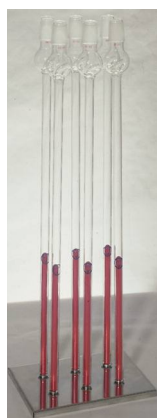

Tube support rack.

Made of chemically treated dural sheet metal, with handles and removable panels that maintain the heat around the tubes and allow the user to examine the samples during the process.

Part No.

- 4005071 for C.O.D. 6 (block 4000507).
- 4005081 for C.O.D. 12 (block 4000508).
- 4005091 for C.O.D. 20 (block 4000509).



Features


Display of time, temperature and program. Selectable temperature from 45 to 450 °C. Memory capacity: 10 programs of 5 stages of each program. Maximum time per stage: 600 minutes. Visual and acoustic alarm at the end of the cycle. Alarm if the temperature probe fails.

Control Panel

- Mains switch.
- Temperature display.
- Push button to select temperature.
- Push button to increase value.
- Push button to start.
- Push button to stop.
- Push button to decrease value.
- Push button to select time.
- Display time and program number
- Push button to select program.

Tubes for C.O.D.

Part No. 1000641

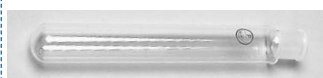

Condenser tubes for C.O.D. process.

Part No. 1000642

